



Formaderm[®]
Dermal Filler Injection



Beauty by

Aid Surg DENTIST
 aAid SurgiAic /ESTAR Gingi
 ivAid Max Forami
 naAid Surg id DENTIST
 Aid SurgiAid DE STAR Gingiv
 Aid SurgiAid DENTIST
 aAid Sur iAid DENTIST R Gingi
 ivAid Yaxi Bone F rami
 naAid irgiAid NTIST
 Aid Surg i DENTIS Gingiv
 Aid Si iAid NTISI
 r GingivAid
 Forma Surgi
 r Gingiv Yaxi Bone
 Max Bone Foramic DENTIST
 r GingivAid
 Forma Surgi
 r GingivAid Max Bone
 r Gingiv Foramic DENTIS
 r GingivAid
 /Ai Yaxi
 imi NTIST
 id rgiA
 3' Bone F
 nic YISIA
 /Aid
 mic 723
 d S
 i
 STA
 /Aid
 axi

INDEX

MAXIGEN BIOTECH INC.

New generation of gentle-style hyaluronate 03
 dermal filler for facial contour restoration

With superior clinical handling property and 04
 long lasting support

Clinical trial and performance 05


Quality system 06






New generation of **gentle-style** hyaluronate dermal filler for contour sculpture

Unique technology of Equilibrium Cross-Linked Hyaluronic Acid

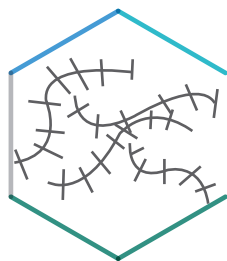

Linear HA

Cross-linking

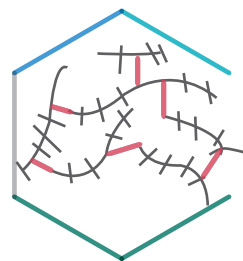


Cross-linked HA

Naturally shaped to create a perfect 3D long lasting support



US FDA audit approved, Taiwan biomedical manufacturer

The expert of hyaluronic acid with mature manufacturing process
Annual sales of over 300 thousand injections
(including Ophthalmic Viscoelastic Device
Intra-articular Injections Dermal Filler and Repair Gel)

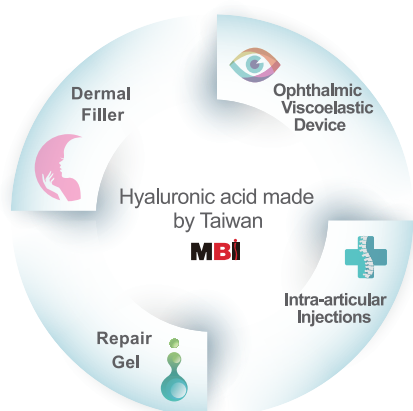


- Unique purification process, eliminating the remain of the BDDE cross-linking agent, high safety and reliability
- Sterilized by high temperature and high pressure procedure, a sterile product with excellent viscoelastic properties

100% purified hyaluronic acid, non animal origin


With superior clinical handling property and providing long lasting support

Superior clinical handling property

Unique technology of Equilibrium Cross-Linked Hyaluronic Acid has created a quality product with perfect cross-linking and uniformity. The smooth and stable extrusion property ensure the precise injections and filling to different parts of facial skin

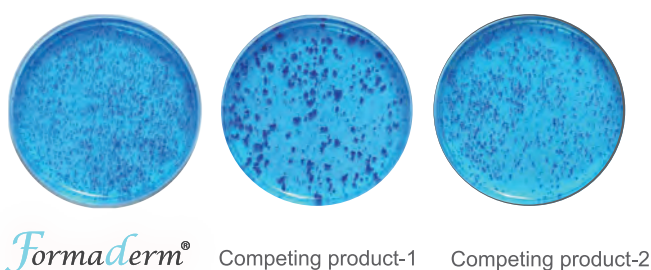

Formaderm® Competing product-1 Competing product-2

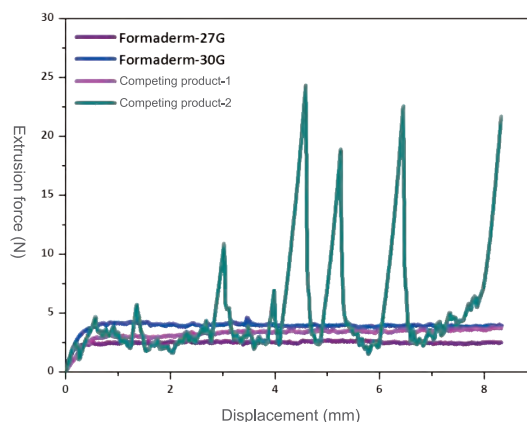

Perfect sustainable support:

Unique technologies of uniformity and cross-linking has created better viscoelastic property than competing products


Samples	G' (Pa) at 5 Hz	G'' (Pa) at 5 Hz
Competing product	445.20 ± 19.78	98.82 ± 7.86
Formaderm®	472.08 ± 33.44	93.04 ± 3.03

Storage Modulus (G') : Elasticity
Loss Modulus (G'') : Viscosity




Clinical trial and performance:


V0 Before injection

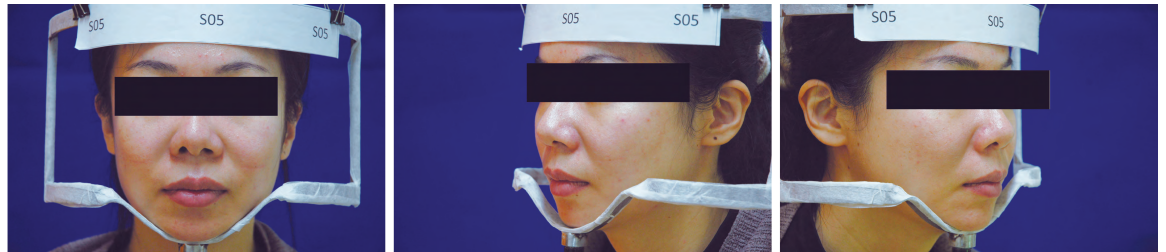

V2 2 weeks after Formaderm(1cc) injection

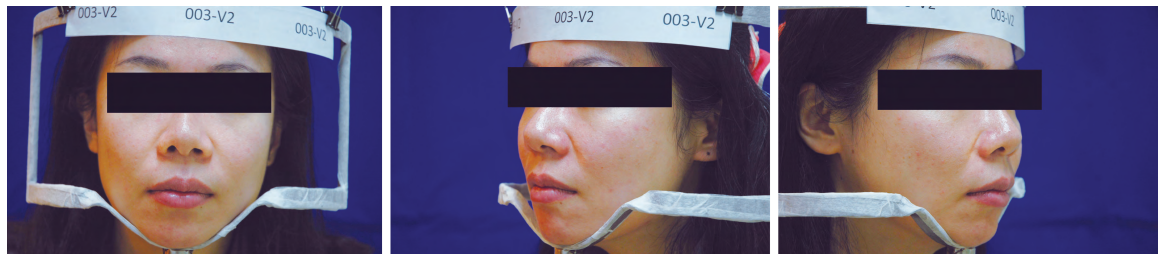

V5 2 4 weeks after Formaderm(1cc) injection

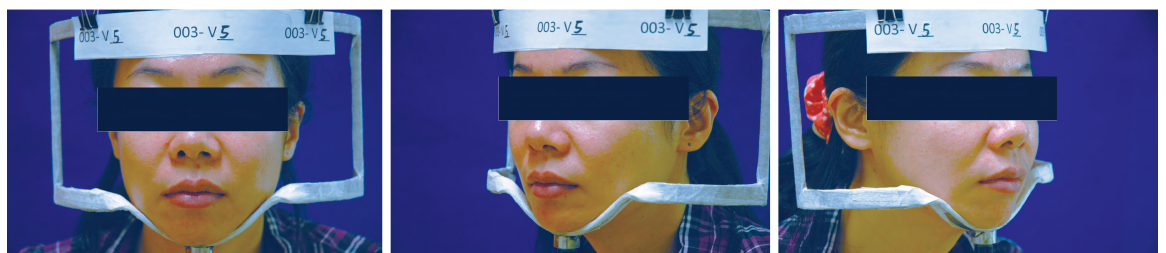

Our clinical results show no significant differences between using Formaderm Dermal Filler and Restylane products in nose augmentation, lip enhancement and wrinkle correction. Both products are effective in facial contour restoration.

Quality system

Quality System (ISO 13485/GMP/21 CFR 820)





Compliant with international standards and regulations of medical devices (CE/FDA/KFDA/HSA)



Formaderm[®]
Dermal Filler Injection

Beauty by



-  With good natural plasticity
-  Creating a perfect 3D long lasting support

MBI MAXIGEN BIOTECH INC.

MAXIGEN BIOTECH INC.
2F, 3F, 4F, 5F, No.88, Keji 1st Rd., Guishan District, Taoyuan City 33383, Taiwan
TEL : +886-3-3287-222 FAX : +886-3-3287-333
www.mbi.com.tw E-MAIL : service@mbi.com.tw

【For medical professional use only】

20170320-D